



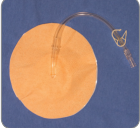




extriCARE™ 2400 Negative Pressure Wound Therapy (NPWT) System Dressings

The extriCARE™ NPWT system incorporates the use of extriCARE™ pump models with Devon Medical Products' proprietary NPWT dressings.

These dressings utilize a wound pad enclosed in mesh under an adhesive film that can be applied directly to the wound bed. A drainage tube connects each dressing with a canister that is connected to the pump. The canister collects exudates and bodily fluid that is removed from the wound bed.

Application of the dressing is very quick and easy, allowing wound care professionals and caregivers to decrease application time.

Dressings are available in various shapes and sizes and designed for use with wounds on all parts of the body.

extriCARE™ Bandages	Type of Dressing	Dressing Size	Gauze Pad Size
	Mini Oval EC2400-MN	120 mm x 130 mm (4.7" x 5.1")	50 mm x 60 mm (2.0" x 2.3")
	Small Oval EC2400-SO	157 mm x 168 mm (6.2" x 6.6")	80 mm x 100 mm (3.1" x 3.9")
	Large Oval EC2400-LO	192 mm x 270 mm (7.6" x 10.6")	77 mm x 152 mm (3.0" x 6.0")
	Sacral Oval EC2400-S	245 mm x 225 mm (9.6" x 8.9")	160 mm x 145 mm (6.3" x 5.7")
	Incision EC2400-IM	140 mm x 370 mm (5.5" x 14.6")	70 mm x 290 mm (2.8" x 11.4")

LB130 Rev A 20121109